

Muscle Sympathetic Nerve Activity Burst Identification and Analysis Standard Operating Procedure

Developed by Miguel Anselmo M.S.

Cardiovascular Research and Rehabilitation Laboratory, University of Minnesota

Table of Contents

Glossary	2
LabChart Setup	2
Cut Files and Prepping Macros	2
Cutting Files	2
Prepping LabChart Cut Files	2
Setting up Excel File for Analysis	10
Performing Burst Identification and Analysis	11
Fixing Errors	16

Glossary

- a. **Continuous Files:** files that do not have segments or intervention periods. For example, a 10-min baseline resting file.
 - b. **CPT:** cold-pressor test
 - c. **ECG:** electrocardiogram
 - d. **EPR:** exercise-pressor reflex
 - e. **MSNA:** muscle sympathetic nerve activity
 - f. **NIBP:** non-invasive blood pressure
 - g. **Segmented Files:** files that have segments or intervention periods. For example, an exercise pressor reflex testing file has four segments (baseline, exercise, post-exercise circulatory occlusion, and recovery).
-

LabChart Setup

IMPORTANT NOTE: It will take time to initially set up this program for analysis, but once all components are created, all subsequent analyses will be very quick. This section will detail, from the beginning, how to analyze continuous baseline files. This program will facilitate the creation of just one workbook for muscle sympathetic nerve activity (MSNA) data per participant to analyze all MSNA data for each study. It is suggested to create a template workbook for each study, so this workbook can be duplicated for each new participant.

Cut Files and Prepping Macros

Setting up LabChart Audio Playback: If this feature is not already installed in LabChart, follow the link below to install Audio Output on LabChart. LabChart may need to be restarted after installation.

<https://www.adinstruments.com/products/audio-output?srsId=AfmBOopGR78B3IYwoWMSAYGGv2qMI2V1P5rFDrUCBgJGzvx-CF-UdrDf>


Cutting Files

- I. Create separate LabChart files for each study procedure of interest. For example, one LabChart file will be baseline rest, another LabChart file will be, for example, a cold-pressor test file.
 - A. For all continuous files (baseline rest, or any files that do not have separate sections such as baseline, stressor, recovery), the file should be highlighted from comment to comment. Next, *File → Save Selection* should be selected, and the selection saved in the desired folder. Reference the ADInstruments website for additional help with this step.
 - B. For all segmented files that have separate sections such as baseline, stressor, recovery, the file should be highlighted from the FIRST comment (for example, Baseline CPT [cold pressor test]) up until the final comment following recovery if that data is of interest (for example, CPT Recovery End). Next, *File → Save Selection* should be selected, and the selection saved in the desired folder. Again, all files should be kept structured and organized.

Prepping LabChart Cut Files

- I. Open the LabChart file that will be analyzed (e.g., SUBID_Rest or SUBID_CPT, etc.).

- II. The time scale should be changed to seconds from the start of the file: **right-click the area at the bottom near the time stamps** → **Time Mode** → **select "Show time from start of file"** → **select "Show time as seconds"**. At the beginning of the file, the time scale should read 0, and the end should read the total? file duration in seconds.
- III. All existing comments should be deleted. If all existing comments are not deleted, downstream errors in the analysis may occur.
IMPORTANT NOTE: If segmented MSNA analysis is being performed (analysis during stressors), the length of each segment (i.e., baseline, CPT, CPT recovery) should be recorded in a separate document **BEFORE** all comments are deleted. This information will be needed in order to appropriately perform the MSNA analysis in the workbook.
- IV. Ensure that the peak of the non-invasive blood pressure (NIBP) waveform is aligned with the electrocardiogram (ECG) peak to ensure that the appropriate diastolic blood pressure (BP) value falls within the correct cardiac cycle. **Click on the dropdown menu shown in the image below** → **Shift** → **shift as necessary depending on the hardware setup**. Also note that different hardware requires different time shifts. For example, if using the reconstructed brachial BP on a Finapres Nova Device, the negative shift is 4s, since the device requires 4s before the reconstructed BP is calculated and sent to LabChart.

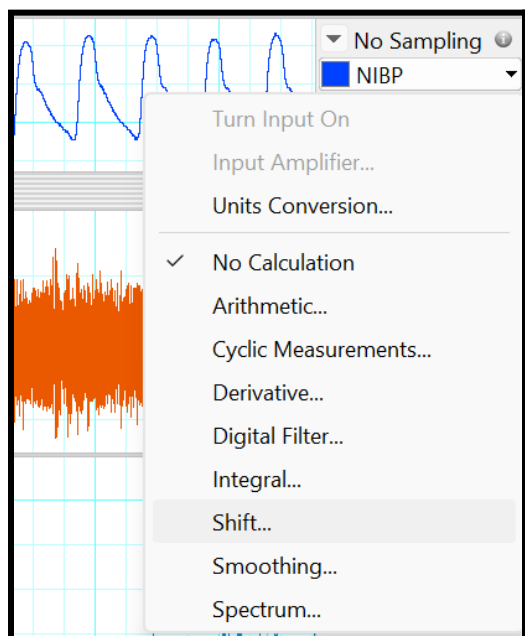

The next steps (VI to XIV) will detail how to set up macros so that these steps do not need to be manually followed with every participant file.

V. Channel Setup:

***IMPORTANT NOTE:** when creating these macros it is suggested to record the steps that are being completed and label the Macro as Channel Setup. If macros are recorded, then they only need to be set up once and they can be uploaded into subsequent files.


To create a macro: **Macro** → **Manage** → **New** → **Name Macro, Show in toolbar** → **Click Record**

ECG: A visible ECG signal is needed for BEAT comments to be placed (next steps).


NIBP: Depending on the hardware, as previously stated, this signal may need to be shifted.

Integrated/RMS MSNA: Either signal will work in this workflow.


Mean Arterial BP: *Cyclic Measurements* → *see image below*

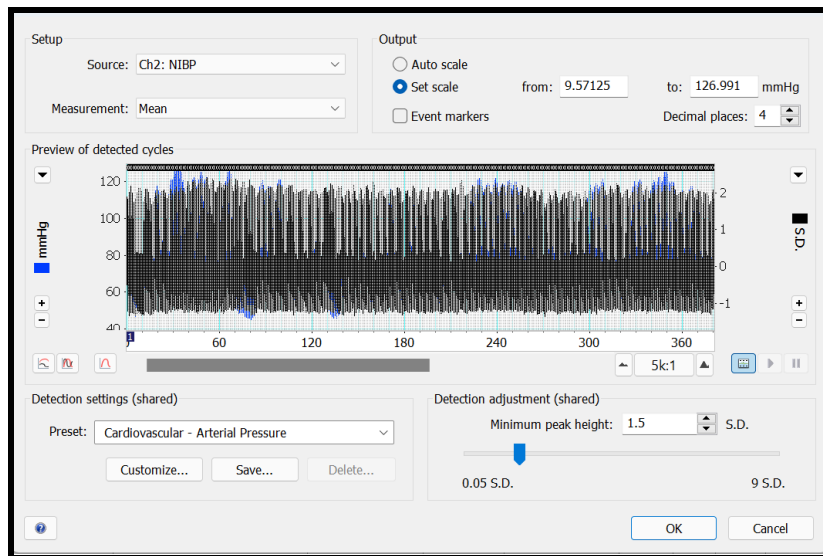

Diastolic BP: *Cyclic Measurements* → see image below

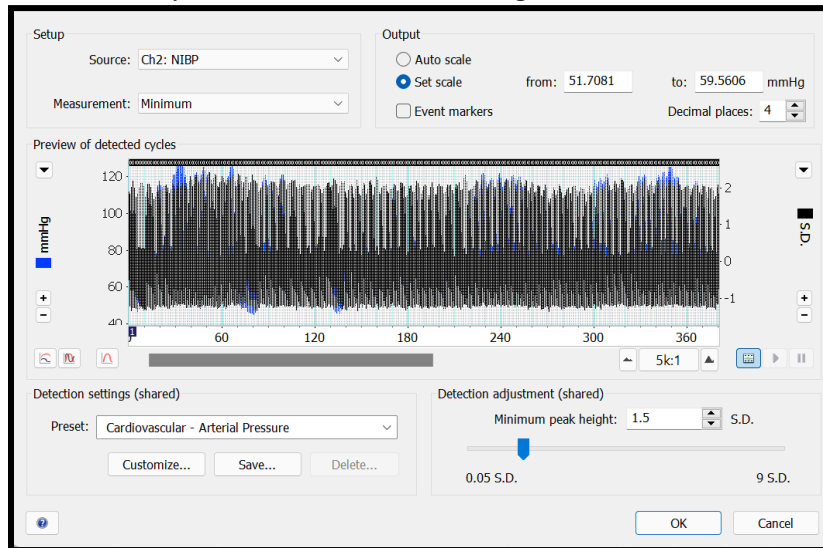

Shifted Diastolic BP: Create a separate channel (right click on screen and press *Add Channel* or go to *Setup* → *Channel Settings* → Increase the # of channels at the bottom of the popup window) → shift → see image below. This channel will be used to calculate sBRS, which incurs a ~1.2s physiological delay, where the burst propagated based on the diastolic BP (DBP) value occurs ~1.2s later. In reality, a range exists here, around 1.1–1.4s, but data exportation is unlikely to be affected because of how the DBP values are extracted. This number may be adjusted as seen fit for a given laboratory. The channel number (Ch6) is dependent on the channel # at which the diastolic BP channel is located.

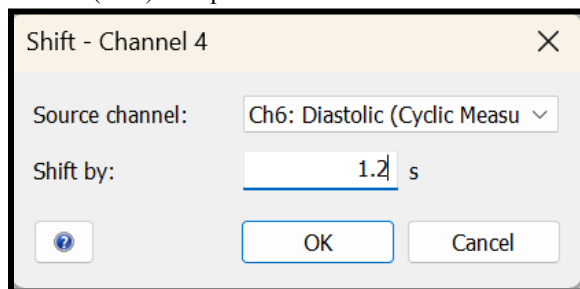

→ **ONCE THESE STEPS ARE COMPLETED PRESS** Stop recording and save.

VI. **Data Pad Setup:** Set up columns in the Data Pad as stated below:

*Tip: record a macro of these steps being completed and label the macro Datapad Setup or something similar.

Macro → **Manage** → **New** → **Name Macro, Show in toolbar** → **Click Record**

**There is no burst area calculation included here, but it is possible. Please refer to the manuscript for limitations regarding burst amplitude and burst area calculations using this program.

Column A: Integrated MSNA Amplitude:

Statistics → Maximum - Minimum → Source Channel: Integrated/RMS MSNA

Column B: Integrated MSNA Comment Number:

Comments → Comment Number → Source Channel: Integrated/RMS MSNA

Column C: Time (s):

Selection & Active Point → Time

Column D: Integrated MSNA Comment Time (s):

Comments → Comment Time → Source Channel: Integrated/RMS MSNA

Column E: ECG Full Comment Text:

Comments → Full Comment Text → Source Channel: ECG

Column F: ECG Comment Number:

Comments → Comment Number → Source Channel: ECG

Column G: Mean Arterial BP Value (mmHg):

Selection & Active Point → Value → Source Channel: Mean Arterial BP

Column H: Diastolic BP Value (mmHg):

Selection & Active Point → Value → Source Channel: Diastolic BP

Column I: Shifted Diastolic BP Value (mmHg):

Selection & Active Point → Value → Source Channel: Shifted Diastolic BP

Column J: Selection Duration (s):

Selection & Active Point → Selection Duration

Turn OFF all other channels by clicking on the column header and selecting Off. Not necessary, but it keeps the Datapad visually clean and easier to export into Excel.

Data Pad headers should look like the image below:

A	B	C	D	E	F	G	H	I	J
Integrated MSNA Maximum - Minimum	Integrated MSNA Comment Number	Time	Integrated MSNA Comment Time	ECG Full Comment Text	ECG Comment Number	Mean Arterial BP Value	Diastolic Value	Shifted Diastolic BP Value	Selection Duration s

→ **ONCE THESE STEPS ARE COMPLETED PRESS Stop recording and save.**

VII. Create Burst ID Macro:

This macro is used to manually mark MSNA bursts.

Macro → Manage → New → Type: Simple Macro, Name: BurstID → Check Show in toolbar → Click OK → Underneath the line Sub Burst ID (), copy the code below:

Call Doc.AddCommentAtInsertionPoint (X, "BURST")

For X, use the n-1 rule to enter in as Channel # for the Integrated/RMS MSNA Channel. For example, if the Integrated/RMS MSNA Channel is Channel 10, then enter the number 9, and the complete code should read:

```

Run Record Stop
Sub BurstID ()
    Call Doc.AddCommentAtInsertionPoint (9, "BURST") ' Creates a comment called 'BURST' in the 10th channel.
End Sub

```

There are automated ways to perform burst marking, which are beyond the scope of the current paper. Please contact the corresponding author if interested in learning about these methods.

VIII. Create a DATAPAD_BEAT Macro:

This macro will look for BEAT comments (these will be placed between each R-R interval in the ECG Channel, and then extract data--the data columns that were set up in the Data Pad-- that are associated with each BEAT comment).

Macro → Manage → New → Type: Simple Macro, Name: DATAPAD_BEAT → Check Show in toolbar → Click OK → Underneath the line Sub DATAPAD_BEAT (), copy the code below:

```

SourceChannelIndex = X
Comment = "beat"
RepeatSetAction = kSelectRegion
RegionSecs = 0.25
SelectMode = kSelectAround
NextComment = ""
SelectScope = kWholeFile
Call Doc.MultipleAddToDataPad_Comment (SourceChannelIndex, Comment, RepeatSetAction, RegionSecs, SelectMode,
NextComment, SelectScope)

```

For X, use the n-1 rule to enter in as the Channel # for the ECG Channel. For example, if the ECG Channel is Channel 1, then enter the number 0, and the complete code should read:

```

Run Record Stop
Sub DATAPAD_BEAT ()
    ' Begin MultipleAddToDataPad_Comment
    SourceChannelIndex = 0
    Comment = "beat"
    RepeatSetAction = kSelectRegion
    RegionSecs = 0.25
    SelectMode = kSelectAround
    NextComment = ""
    SelectScope = kWholeFile
    Call Doc.MultipleAddToDataPad_Comment (SourceChannelIndex, Comment, RepeatSetAction, RegionSecs, SelectMode, NextComment, SelectScope)
    ' End MultipleAddToDataPad_Comment
End Sub

```

IX. Create a DATAPAD_BURST Macro:

This macro will look for BURST comments (the comments made when bursts are manually/automatically marked in the Integrated/RMS MSNA channel). Data (i.e., the data columns set up in the Data Pad) associated with each BURST comment will then be extracted. In LabChart, data needs to be exported based on BURST comments separately from BEAT comments, because all bursts in a file cannot be time-aligned with BEAT comments. Thus, in order to preserve data associated with burst amplitude, these data need to be exported separately. This macro will pull data from a 0.5s. search window surrounding each BURST comment.


Macro → Manage → New → Type: Simple Macro, Name: DATAPAD_BURST → Check Show in toolbar → Click OK → Underneath the line Sub DATAPAD_BURST (), copy the code below:

```

SourceChannelIndex = X
Comment = "BURST"
RepeatSetAction = kSelectRegion
RegionSecs = 0.5
SelectMode = kSelectAround
NextComment = ""
SelectScope = kWholeFile
Call Doc.MultipleAddToDataPad_Comment (SourceChannelIndex, Comment, RepeatSetAction, RegionSecs, SelectMode,
NextComment, SelectScope)

```

For **X**, use the *n-1* rule to enter in as the Channel # for the Integrated/RMS MSNA Channel. For example, if the Integrated/RMS MSNA Channel is Channel 11, then enter the number 10, and the complete code should read:

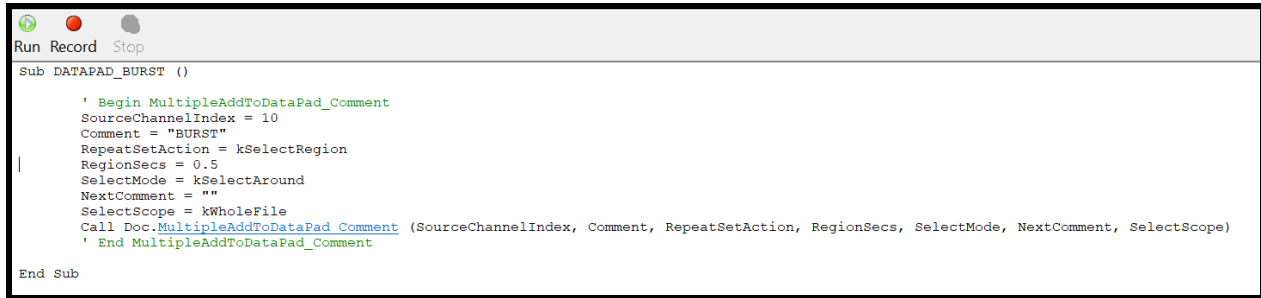

```

Sub DATAPAD_BURST ()

    ' Begin MultipleAddToDataPad_Comment
    SourceChannelIndex = 10
    Comment = "BURST"
    RepeatSetAction = kSelectRegion
    RegionSecs = 0.5
    SelectMode = kSelectAround
    NextComment = ""
    SelectScope = kWholeFile
    Call Doc.MultipleAddToDataPad_Comment (SourceChannelIndex, Comment, RepeatSetAction, RegionSecs, SelectMode, NextComment, SelectScope)
    ' End MultipleAddToDataPad_Comment

End Sub

```

X. Create a Macro automating BEAT comments in the ECG Channel

This step requires the creation of several macros: (1) ECG_MAKESEL & (2) MARK_ECG

This step will create the **ECG_MAKESEL** Macro

Macro → Manage → New → Type: Simple Macro, Name: ECG_MAKESEL → Show in toolbar → Click OK → Underneath the line Sub ECG_MAKESEL (), copy the code below:

```

ChannelIndex = X
SetAction = kSetEndOfSelection
SelectMode = kSelectAround
SelectTime = 1
DataDisplayMode = kViewDataVisible
SelectAll = False
Direction = kSearchForward
FindType = "Move backward"
FindData = "AmountToMove=0.1;"
Call Doc.Find (ChannelIndex, SetAction, SelectMode, SelectTime, DataDisplayMode, SelectAll, Direction, FindType, FindData)

```

```

ChannelIndex = X
SetAction = kSetActivePoint
SelectMode = kSelectAround
SelectTime = 1
DataDisplayMode = kViewDataVisible
SelectAll = False
Direction = kSearchForward
FindType = "Local maxima"
FindData = "NoiseThreshold=0.05;"
Call Doc.Find (ChannelIndex, SetAction, SelectMode, SelectTime, DataDisplayMode, SelectAll, Direction, FindType, FindData)

```

```


ChannelIndex = X
SetAction = kSetToPreviousPoint
SelectMode = kSelectAround
SelectTime = 1
DataDisplayMode = kViewDataVisible
SelectAll = False
Direction = kSearchForward
FindType = "Local maxima"
FindData = "NoiseThreshold=0.05;"
Call Doc.Find (ChannelIndex, SetAction, SelectMode, SelectTime, DataDisplayMode, SelectAll, Direction, FindType, FindData)

```

The NoiseThreshold is set to 0.05 in this example. Depending on the quality and/or presence of artifacts in the ECG signal, the threshold may need to be adjusted. More detail is provided in a subsequent section, "Performing Burst Identification and

Analysis."

For **X**, use the $n-1$ rule to enter in as the Channel # for the ECG Channel. For example, if the ECG Channel is Channel 1, then enter the number 0, and the complete code should read:

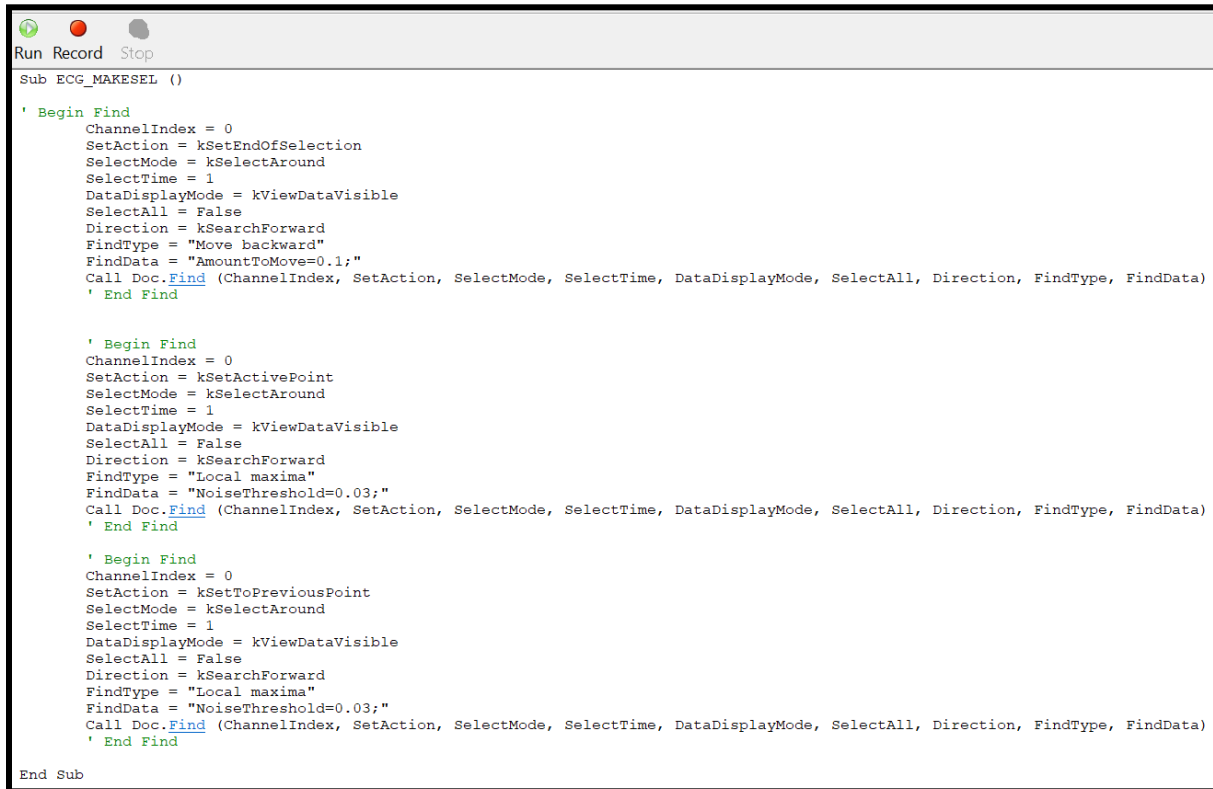

```
Run Record Stop
Sub ECG_MAKESEL ()

' Begin Find
ChannelIndex = 0
SetAction = kSetEndOfSelection
SelectMode = kSelectAround
SelectTime = 1
DataDisplayMode = kViewDataVisible
SelectAll = False
Direction = kSearchForward
FindType = "Move backward"
FindData = "AmountToMove=0.1;"
Call Doc.Find (ChannelIndex, SetAction, SelectMode, SelectTime, DataDisplayMode, SelectAll, Direction, FindType, FindData)
' End Find

' Begin Find
ChannelIndex = 0
SetAction = kSetActivePoint
SelectMode = kSelectAround
SelectTime = 1
DataDisplayMode = kViewDataVisible
SelectAll = False
Direction = kSearchForward
FindType = "Local maxima"
FindData = "NoiseThreshold=0.03;"
Call Doc.Find (ChannelIndex, SetAction, SelectMode, SelectTime, DataDisplayMode, SelectAll, Direction, FindType, FindData)
' End Find

' Begin Find
ChannelIndex = 0
SetAction = kSetToPreviousPoint
SelectMode = kSelectAround
SelectTime = 1
DataDisplayMode = kViewDataVisible
SelectAll = False
Direction = kSearchForward
FindType = "Local maxima"
FindData = "NoiseThreshold=0.03;"
Call Doc.Find (ChannelIndex, SetAction, SelectMode, SelectTime, DataDisplayMode, SelectAll, Direction, FindType, FindData)
' End Find

End Sub
```

XI. This step will create the next macro needed to mark BEAT comments in the ECG channel: MARK_ECG.

Macro → Manage → New → Type: Simple Macro, Name: MARK_ECG → Show in toolbar → Click OK → Underneath the line Sub MARK_ECG (), copy the code below:

```
ChannelIndex = X
SetAction = kSetActivePoint
SelectMode = kSelectAround
SelectTime = 1
DataDisplayMode = kViewDataVisible
SelectAll = True
Direction = kSearchForward
FindType = "Start of file"
FindData = ""
Call Doc.Find (ChannelIndex, SetAction, SelectMode, SelectTime, DataDisplayMode, SelectAll, Direction, FindType, FindData)
startRec = Doc.SelectionStartRecord
Do While (IsCurrentSelectionWithin (startRec, 0, startRec, kRecordEndOffset) = True)
Call Doc.Macros.ECG_MAKESEL
If (Services.ShouldExitCurrentRepeat()) Then Exit Do
If Not (IsCurrentSelectionWithin (startRec, 10, startRec, kRecordEndOffset) = True) Then Exit Do
Call Doc.AddCommentAtInsertionPoint (X, "BEAT")
If Not (IsCurrentSelectionWithin (startRec, 0, startRec, kRecordEndOffset) = True) Then Exit Do
Loop
If Not (startRec = kEndRecordIndex) Then
Call Doc.SetSelectionRange (startRec, 2147483647, startRec, 2147483647)
End If

End Sub
```


This macro “calls” the **ECG_MAKESEL** macro, and essentially uses this macro to find local maxima (or each R-wave in the ECG tracing), and then places a BEAT comment between each R-wave.

For **X**, use the *n-1* rule to enter in as the Channel # for the ECG Channel. For example, if the ECG Channel is Channel 1, then enter the number 0, and the complete code should read:

```

Run Record Stop
Sub MARK_ECG ()

    ' This example runs through a block, recording local maxima.
    '

    ChannelIndex = 0
    SetAction = kSetActivePoint
    SelectMode = kSelectAround
    SelectTime = 1
    DataDisplayMode = kViewDataVisible
    SelectAll = True
    Direction = kSearchForward
    FindType = "Start of file"
    FindData = ""
    Call Doc.Find (ChannelIndex, SetAction, SelectMode, SelectTime, DataDisplayMode, SelectAll, Direction, FindType, FindData)

    startRec = Doc.SelectionStartRecord

    Do While (IsCurrentSelectionWithin (startRec, 0, startRec, kRecordEndOffset) = True)
        Call Doc.Macros.ECG_MAKESEL

        If (Services.ShouldExitCurrentRepeat()) Then Exit Do
        If Not (IsCurrentSelectionWithin (startRec, 10, startRec, kRecordEndOffset) = True) Then Exit Do


        Call Doc.AddCommentAtInsertionPoint (0, "BEAT")
        If Not (IsCurrentSelectionWithin (startRec, 0, startRec, kRecordEndOffset) = True) Then Exit Do
    Loop
    If Not (startRec = kEndRecordIndex) Then
        Call Doc.SetSelectionRange (startRec, 2147483647, startRec, 2147483647)
    End If

End Sub


```

XII. **Record a macro to automate setting up the channel view.** The key channels are the ECG, NIBP, raw or filtered MSNA (depending on preference), and the integrated/RMS MSNA channel.

Macro → Manage → New → Type: Simple Macro, Name: Channel View → Show in toolbar → Click OK → Record setting up the channel view, similar to the image below:

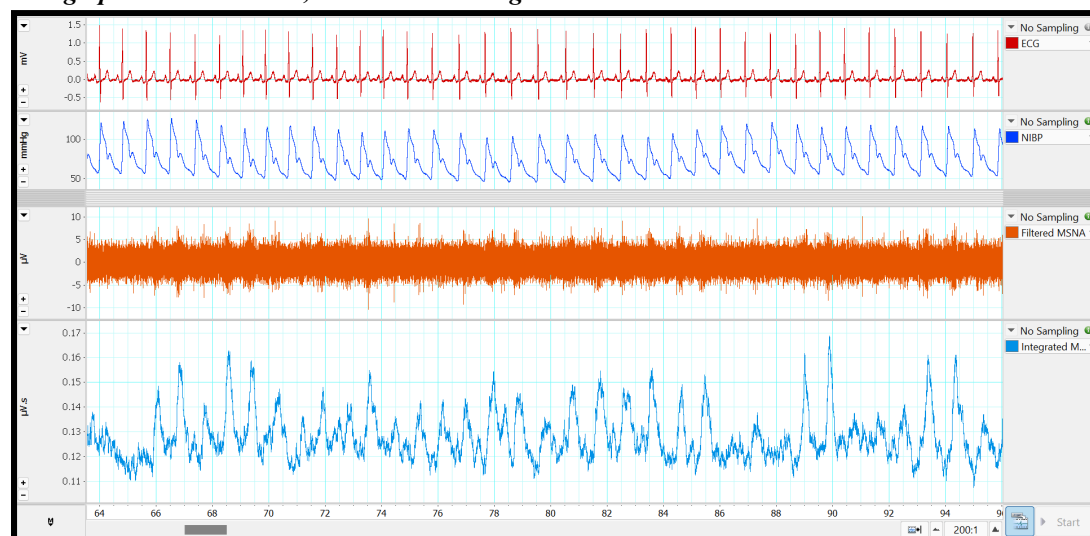

The setup macro list should look similar to this:

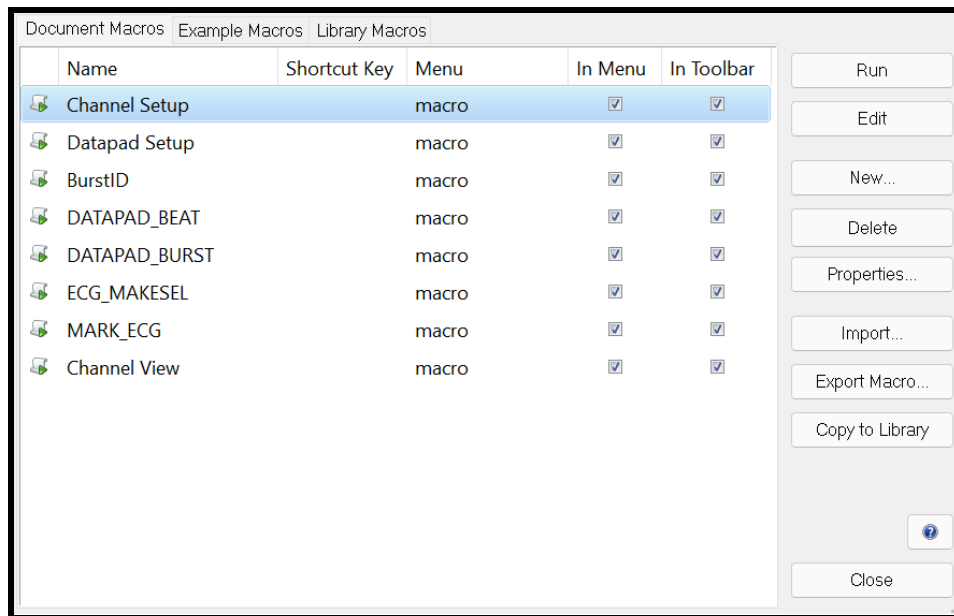

- XIII. **TIP!:** Combine macros to expedite the setup process. Create a new macro, *run* the macros that will be combined, save as a parent macro.

Setting up Excel File for Analysis

- I. Open the BurstID_TemplateFile on Excel. Make sure Excel is allowed to run macros. A popup in Excel may occur if macros need to be enabled for this file.
- II. In the CumulativeSheet tab, enter the specified participant's SUBID and Group # (if applicable; for example, females are Group 1 and males are Group 2):

	A	B
1	SUBID	Group
2	Example001	1
3		

- III. Excel Workbook Tabs:

CumulativeSheet: all data will automatically populate into here from subsequent sheets

LabChart_CorrectedInput_Rest: this sheet will calculate resting MSNA variables

LabChartInput_BURST_Rest: data here will be pasted from the Data Pad after extraction based on BURST comments from resting LabChart files

LabChartInput_BEAT_Rest: data here will be pasted from the Data Pad after extraction based on BEAT comments from resting LabChart files

MSNA Outcomes_Rest: all resting MSNA outcomes will be populated here for ease of viewing

BurstTransduction_MAP: resting MAP transduction data will be populated here

NonBurstTransduction_MAP: resting MAP nonburst transduction data will be populated here

BurstTransduction_DBP: resting DBP transduction data will be populated here

NonBurstTransduction_DBP: resting DBP nonburst transduction data will be populated here

CPT/Other Stressor Tabs: While these tabs are set up to analyze a stressor with a 3-min baseline, 2-min stressor, and 5-min recovery, these tabs are modifiable to accommodate different durations of analyses. If clarity is needed on how to adjust formulas and outcome variable nomenclature, please contact the corresponding author.

MSNA Outcomes_CPT: all MSNA outcomes from CPT or other stressors will be populated here for ease of viewing. Modify this sheet as needed to analyze various stressors (Stroop Color and Word Test, mental arithmetic, etc.) that contain a baseline, stressor, and recovery.

LabChart_CorrectedInput_CPT: this sheet will calculate MSNA variables for CPT/other stressors

LabChartInput_BURST_CPT: data here will be extracted based on BURST comments from CPT/other stressor LabChart files

LabChartInput_BEAT_CPT: data here will be extracted based on BEAT comments from CPT/other stressor LabChart files

EPR (Exercise Pressor Reflex) Stressor Tabs:

While these tabs are set up to analyze exercise pressor reflex testing (3-min baseline, 2-min handgrip [IHG], 3-min post-exercise circulatory occlusion [PECO]), this tab is modifiable to accommodate different durations of analyses. For example, the IHG outcomes are calculated as 25% of the IHG duration, and can thus accommodate a wide range of exercise times. However, PECO is calculated in 30-ss epochs, but can accommodate a wide range of durations. If clarity is needed on how to adjust formulas and outcome variable nomenclature, please contact the corresponding author.

MSNA Outcomes_EPR: all exercise pressor reflex testing MSNA outcomes will be populated here for ease of viewing. This sheet can be used for any stressor with four segmented parts. Modify this sheet as needed.

LabChart_CorrectedInput_EPR: this sheet will calculate EPR/other stressor MSNA variables


LabChartInput_BURST_EPR: data here will be pasted from the Data Pad after extraction based on BURST comments from EPR/other stressor LabChart files

LabChartInput_BEAT_EPR: data here will be pasted from the Data Pad after extraction based on BEAT comments from EPR/other stressor LabChart files

IV. Move onto the next section.

Performing Burst Identification and Analysis

- I. **IMPORTANT NOTE:** Save the current file with all the created macros so that they can be reuploaded into all subsequent files. In new files choose: *Macro* → *Manage* → *Import* → *Navigate to saved macro file and import into new file*.
- II. Run the **Channel Setup** macro.
- III. Run the **Datapad Setup** macro.
- IV. Click anywhere on the ECG tracing. Run the **ECG_MAKESEL** macro a few times. It should highlight from each R-R interval on the ECG tracing as seen below:

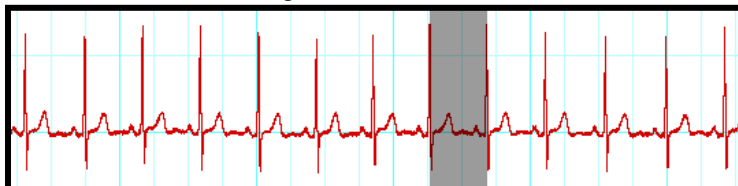

If this macro is skipping beats, or erroneously selecting other components of the tracing that are not the R-wave, stop the

macro, delete all of the BEAT comments, and troubleshoot by editing certain lines of the code:


Option 1: Either increase or decrease the noise threshold. Be strategic, and notice what increasing and decreasing this number does to detecting components of the tracing. Remember, there are two places in the ECG_MAKESEL macro where this number will need to be changed.

Option 2: Change the macro to detect the local minimas rather than local maximas.


FindType = "Local maxima" change to FindType = "Local minima"

If neither of these options work for all beats, it may be necessary to run the MARK_ECG macro, then manually enter some BEAT comments or to move the placement of some BEAT comments. Once all R-R intervals have been marked by a BEAT comment, proceed to the next step. During data collection, ensuring that there is a clean ECG signal will reduce the likelihood of having to manually enter and/or edit BEAT comments placed by the macro.


- V. Double-click on the timestamp scale at the bottom of the screen to highlight the entire file:

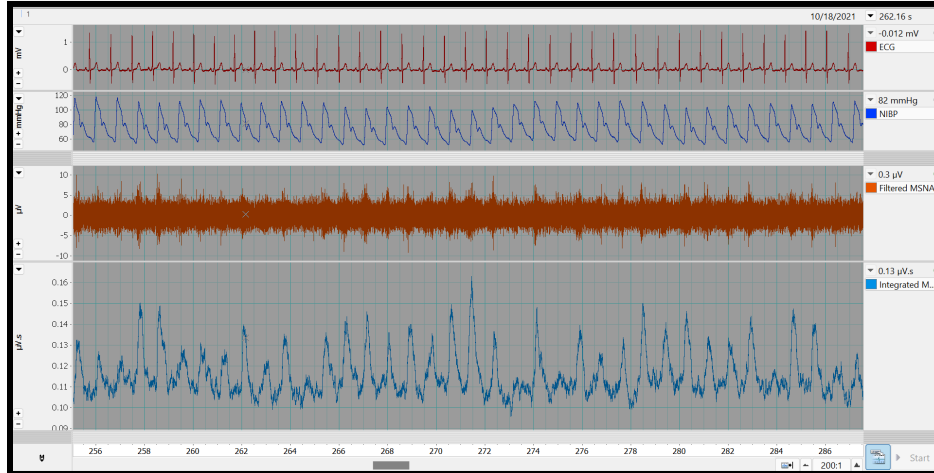

- VI. Run the MARK_ECG Macro. This macro will add a BEAT comment between each R-wave:

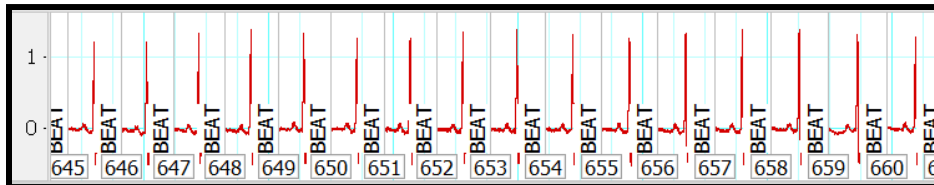

- VII. It is now time to perform burst identification. In the integrated MSNA/RMS MSNA, click on the peak of each burst, then run the Manual_BurstID macro:

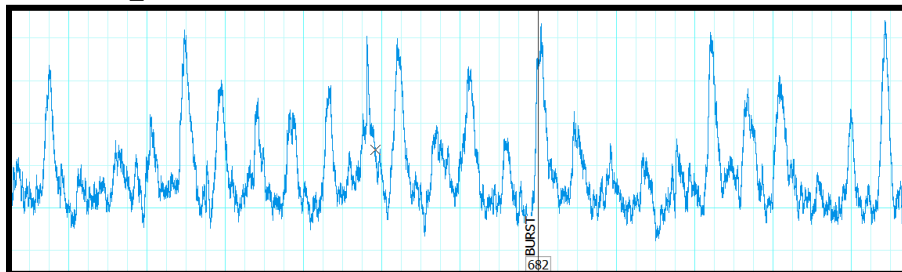

This will add a BURST comment at the peak of each burst. **IMPORTANT NOTE:** (1) automated methods may be developed so that each burst does not need to be manually annotated, and (2) the duration from the preceding R-wave may be checked to ensure that bursts fall within the normal range of burst propagation. The numerical order of the BURST comments will not affect subsequent data analysis, so comments do not need to be shifted if an additional BURST comment is added between two existing BURST comments. The steps detailed below may be of assistance toward burst identification:

1. Assess file quality. Determine whether the file is clean or messy (i.e bursts that are audible with correct morphology that may be shrouded by noise, skin sympathetic activity present, etc.).

2. Evaluate burst morphology and signal-to-noise ratio. A signal-to-noise ratio below 3:1 does not, on its own, disqualify a burst, morphology and the criteria below should be weighed together.

3. Resolve uncertain bursts. If a potential burst may be obscured by electrical noise, use audio playback to listen for it and check the R-wave latency. Bursts typically fall 1-1.4s after the cardiac cycle that generated them.

4. Cross-reference channels. Use the raw and filtered MSNA channels to gain additional information about bursts appearing in the integrated/RMS MSNA channel. If there is a burst present, sloped upward and downwards deflections typically appear in the filtered signal when there is a burst.

Once burst marking for the entire file is completed, move onto the next step.

VIII. Open the **BurstID_TemplateFile_Date*** Excel workbook. Allow macros to be run in the Excel Workbook after downloading and opening the file. There should be no security risks to do this.

IX. Go to the **LabChartInput_BURST_Rest** tab.

X. Run the **DATAPAD_BURST** macro. This macro will look for BURST comments and add all pertinent, associated data into the Data Pad. Using the “Cell edit mode,” copy columns A-I of the Data Pad (not the column names) into the **LabChartInput_BURST_Rest** tab:

Step 1: Add all data into Data Pad from BURST comments by running the DATAPAD_BURST macro.

	A Integrated MSNA Max - Min µV.s	B Integrated MSNA Comment Number	C Time s	D Integrated MSNA Cmt Time s	E ECG Full Comment Text	F ECG Comment Number	G Mean Arterial BP Value mmHg	H Diastolic Value mmHg	I Shifted Diastolic BP Value mmHg	J Selection Duration s
1	0.021794193	1224	599.3057	599.4307	BEAT	680	82.01477051	57.88750076	57.92499917	0.25
2	Integrated MS	Integrated MS		Integrated MS	ECG	ECG	Mean Arterial	Diastolic	Shifted Diast	
3	Max - Min	Cmt No	Time	Cmt Time	Cmt Text	Cmt No	Value	Value	Value	Sel Duration
4	µV.s		s	s			mmHg	mmHg	mmHg	s
5	0.034667356	682	0.2197	0.3447			#NUM!	#NUM!	#NUM!	0.25
6	0.029755910	683	0.845	0.97			#NUM!	#NUM!	#NUM!	0.25
7	0.024505513	684	2.6362	2.7612	BEAT	2	72.56070709	49.91249847	50.91249845	0.25
8	0.016746868	685	3.6829	3.8079	BEAT	3	71.85609436	48.91249847	49.91249843	0.25
9	0.022764071	686	4.7093	4.8343	BEAT	4	69.59938812	48.20624924	48.9124984	0.25
10	0.027770381	687	5.527	5.652	BEAT	5	66.10929108	46.89374924	48.20624922	0.25
11	0.015461599	688	6.5623	6.6873	BEAT	6	66.81414032	47.91249847	46.89374945	0.25
12	0.014867998	689	7.5419	7.6669	BEAT	7	69.11362457	48.91249847	47.91249838	0.25
13	0.012768098	690	8.4861	8.6111	BEAT	8	70.72471619	48.85625076	48.9124984	0.25
14	0.032080768	691	9.4042	9.5292	BEAT	9	70.30609894	48.19374847	48.85625074	0.25
15	0.037842417	692	10.3375	10.4625	BEAT	10	69.44993591	48.88750076	48.19374867	0.25
16	0.035301492	693	11.3095	11.4345	BEAT	11	66.83261108	48.89374924	48.88750067	0.25
17	0.014983729	694	12.3731	12.4981	BEAT	12	70.06891632	50.93125153	48.89374902	0.25
18	0.047340212	695	13.4955	13.6205	BEAT	14	74.58034515	52.17499924	50.93125167	0.25
19	0.020510838	696	14.2387	14.3637	BEAT	14	74.58034515	52.17499924	51.88750074	0.25

Step 2: Copy data AS SHOWN ABOVE into the LabCharInput_BURST_Rest tab

Integrated MSNA Amplitude	Integrated MSNA Comment #	Time (s)	Integrated MSNA Comment Time (s)	ECG Comment Text	ECG Comment #	Mean Arterial BP Value mmHg	Diastolic BP Value mmHg	Shifted Diastolic BP Value mmHg
0.034667267	682	0.22	0.34			#NUM!	#NUM!	#NUM!
0.02975591	683	0.85	0.97			#NUM!	#NUM!	#NUM!
0.024505514	684	2.64	2.76 BEAT		2	72.56070709	49.91249847	50.91249845
0.016746869	685	3.68	3.81 BEAT		3	71.85609436	48.91249847	49.91249843
0.022764071	686	4.71	4.83 BEAT		4	69.59308812	48.20624924	48.9124984
0.027770382	687	5.53	5.65 BEAT		5	66.10929108	46.89374924	48.20624922
0.0154616	688	6.56	6.68 BEAT		6	66.81414032	47.91249847	46.89374945
0.014867998	689	7.54	7.67 BEAT		7	69.10362457	48.91249847	47.91249838
0.012768098	690	8.49	8.61 BEAT		8	70.72471619	48.85625076	48.9124984
0.020280768	691	9.40	9.53 BEAT		9	70.30609894	48.19374947	48.85625074
0.037842417	692	10.34	10.46 BEAT		10	69.44992591	48.08750076	48.19374967
0.035301493	693	11.31	11.43 BEAT		11	66.8326108	48.89374924	48.88750067
0.014963729	694	12.37	12.50 BEAT		12	70.0689632	50.9325153	48.89374902
0.047340213	695	13.50	13.62 BEAT		14	74.58034515	52.17499924	50.9325167
0.020510038	696	14.24	14.36 BEAT		14	74.58034515	52.17499924	51.88750074
0.026837458	697	15.22	15.34 BEAT		15	73.88000574	51.91249847	52.17499915
0.02121877	698	16.15	16.27 BEAT		16	73.27601624	52.19374947	51.91249848
0.020462424	699	17.13	17.25 BEAT		18	71.2552948	52.90625	52.19374853
0.020306337	700	18.04	18.17 BEAT		19	74.90862274	54.88750076	52.90624991
0.017680418	701	19.94	20.06 BEAT		21	80.05818939	55.9325153	56.91249848
0.013161248	702	20.74	20.86 BEAT		21	80.05818939	55.9325153	56.93125158
0.01891146	703	21.71	21.83 BEAT		23	78.16384125	55.91249847	55.93125152
0.037588079	704	22.60	22.72 BEAT		24	76.05962971	54.91249847	55.91249849
0.017654515	705	23.42	23.55 BEAT		24	75.43449402	53.21250153	54.91249849
0.02305006	706	25.26	25.38 BEAT		26	75.46483612	51.91249847	53.87499998
0.036359621	707	26.16	26.29 BEAT		27	74.02708435	52.24375153	51.91249848

Step 3: There may be #NUM cells that are highlighted in yellow. This occurs due to the way LabChart clips exported data at the beginning and end of files. Press the button to the right to clear all of these cells.

Integrated MSNA Amplitude	Integrated MSNA Comment #	Time (s)	Integrated MSNA Comment Time (s)	ECG Comment Text	ECG Comment #	Mean Arterial BP Value mmHg	Diastolic BP Value mmHg	Shifted Diastolic BP Value mmHg
0.034667267	682	0.22	0.34					
0.02975591	683	0.85	0.97					
0.024505514	684	2.64	2.76 BEAT					
0.016746869	685	3.68	3.81 BEAT					
0.022764071	686	4.71	4.83 BEAT					
0.027770382	687	5.53	5.65 BEAT					
0.0154616	688	6.56	6.68 BEAT					
0.014867998	689	7.54	7.67 BEAT					
0.012768098	690	8.49	8.61 BEAT					
0.020280768	691	9.40	9.53 BEAT					
0.037842417	692	10.34	10.46 BEAT					
0.035301493	693	11.31	11.43 BEAT					
0.014963729	694	12.37	12.50 BEAT					
0.047340213	695	13.50	13.62 BEAT					
0.020510038	696	14.24	14.36 BEAT					
0.026837458	697	15.22	15.34 BEAT					
0.02121877	698	16.15	16.27 BEAT					
0.020462424	699	17.13	17.25 BEAT					
0.020306337	700	18.04	18.17 BEAT					
0.017680418	701	19.94	20.06 BEAT					
0.013161248	702	20.74	20.86 BEAT					
0.01891146	703	21.71	21.83 BEAT					
0.037588079	704	22.60	22.72 BEAT					
0.017654515	705	23.42	23.55 BEAT					
0.02305006	706	25.26	25.38 BEAT					
0.036359621	707	26.16	26.29 BEAT					

- XI. **COMPLETELY** clear the Data Pad in LabChart. Then, highlight the entire file, and run the **DATAPAD_BEAT** macro. This macro will look for BURST comments and add all pertinent, associated data into the Data Pad. Copy the Data Pad (not the column names) into the **LabChartInput_BEAT_Rest** tab. Press the button to eliminate all #NUM! values.
- XII. Now go to the **LabChart_CorrectedInput** tab. This tab will calculate many resting MSNA outcomes. Scan the tab for any error cells and clear or modify, if applicable. Go to the **MSNA_OutcomesRest** tab. If any data are blank, then that means there may be some errors in either the way the data was exported into LabChart or in the **LabChart_CorrectedInput_Rest** sheet. Please double check that the data is being calculated as intended and that data were exported into the correct columns in the **LabChartInput_BEAT_Rest** and **LabChartInput_BURST_Rest** tabs. The Unweighted and Weighted sBRS slopes will only appear if their respective correlation coefficients are ≥ 0.7 . An example of this tab is shown below:

	A	B
1		
2	MSNA Outcome Variables	
3	Burst Incidence (bursts/100Hb)	72.98091043
4	Burst Frequency (bursts/min)	49.75971165
5	Burst Amplitude (a.u.)	3.289091416
6	Total Activity (a.u.)	163.6642405
7	Unweighed Sympathetic BRS (bursts/100Hb/mmHg)	-3.553133073
8	Weighed Sympathetic BRS (bursts/100Hb/mmHg)	-2.874724816
9		

- XIII. Go to the Transduction tabs. Check for any errors or missing data and modify, if applicable. These tabs will calculate overall transduction, overall nonburst transduction, burst amplitude quartile transduction, and BP tertile transduction. **If there were any baseline shifts or movement of the microelectrode, burst amplitude analysis may not be accurate.**
- XIV. If using the CPT or EPR (other stressor) tabs, follow the process outlined above for the specified BEAT or BURST data. Ensure, in their respective CorrectedInput tabs, to input the specified timepoints/information in the cells that are colored **green**. Use these tabs to perform analyses on a wide range of stressors. Be familiar with the structure and some of the code before making any changes.

N	O	P	Q	R	S	T
	MSNA OUTCOMES CALCULATIONS					
	BL CPT START TIME	0				
	CPT START TIME			<--- Fill S3 out with the time (in seconds) that CPT started (use timestamp in LabChart file).		
	CPT END TIME					
	CPT DURATION			<--- Fill S5 out with the time (in seconds) that CPT lasted (CPT duration).		
	BL CPT Average DBP					
	BL CPT Average MAP					
	CPT Average DBP					
	CPT Average MAP					
	CPT Burst Count					
	CPT Recording Time					
	CPT Total Cardiac Cycles					
	CPT Average BIN					
	CPT Average BF					
	CPT Average Burst Amplitude					
	CPT Average TA					

N	O	P	Q	R
MSNA OUTCOMES CALCULATIONS				
BL IHG START TIME		0		
IHG START TIME				<--- Fill out S3 with the time (in seconds) that IHG started.
IHG END TIME				
IHG TTF				<--- Fill out S5 with the time (in seconds) that that IHG lasted (time to task failure).
PECO START TIME				
PECO END TIME				
PECO DURATION				<--- Fill out S8 with the time (in seconds) that PECO lasted.
BL IHG Average DBP				
BL IHG Average MAP				
IHG Average DBP				
IHG Average MAP				
IHG Burst Count				
IHG Recording Time				
IHG Total Cardiac Cycles				
IHG Average BIN				
IHG Average BF				
IHG Average Burst Amplitude				
IHG Average TA				

Fixing Errors

Follow these steps if any errors arise in completing this analysis:

1. Ensure that data is copied appropriately from LabChart. Make sure all BEAT comments are **only** attached to the ECG channel and all BURST comments are **only** attached to the Integrated/RMS MSNA Channel. Ensure that all pre-existing comments are deleted from LabChart.
2. Make sure the timestamp is changed appropriately, as previously described, to be in seconds from the start of the file.
3. Ensure that the Data Pad columns are accurately aligned with the Excel columns.
4. Some of the formulas in certain cells in the workbook are conditional. Blank cells should be checked to ensure the reference cells are not blank themselves. For example, in the image below, the data present in the DBP Floor (Cell P2) and DBP Ceiling (Cell P3) cell formulas are conditional on Cell J10 having a value. Thus, if Cell J10 is blank (depending on the dataset), this condition should either be deleted, or another cell that has a value that was exported from LabChart should be selected (i.e., J10 should be changed to the number of a different cell that has a value). In the MSNA_Outcomes tabs for rest and during stressors, many of the calculated cells to display analyzed MSNA data have an IFERROR function, so if the cells are blank it means there is an error in the cells that are being

calculated.

=IF(J10="", "", FLOOR(MIN(H5:H1001), 3))				
L	M	N	O	P
Calculations				
DBP Floor	0			
DBP Ceiling				1
Average DBP	0			2
	0			3

5. For the sympathetic baroreflex calculation, incorrect pressure bins may be generated if there are erroneous data present in Column H, the shifted DBP column (see image below; a DBP of -6 mmHg is not possible). If such bins occur, the button titled "Find DBP Values Associated with Erroneous DBP Bins!" should be clicked and the Bin # entered. In Column H (see figure below), the DBP values associated with each bin will be highlighted. If any erroneous DBP values are noticed (such as in the image below, DBP cannot physiologically be a negative value), the highlighted cells should be deleted or troubleshooted further. This will correct the sBRS calculation. Additionally, while the midpoint DBP of each bin is presented in this calculation, in reality, using the floor, midpoint, or ceiling of each bin for each X-value does not matter because it will not influence the slope or the correlation coefficient if they are all shifted equally. The second figure below shows what happens when the macro button is pressed and the bin number is selected (Bin #1 in this example). Easily find in Column H the erroneous blood pressure value and decide whether or not to delete this value.

SYMPATHETIC BAROREFLEX SENSITIVITY CALCULATIONS							
DBP Floor	-6	DBP Bin #	DBP Bins	# Beats/Bin	# Bursts/Bin	Burst Incidence/Bin	sBRS (unweighted)
DBP Ceiling	72	#1	-6	1	1	100	-0.517058214
Average DBP	58.76146705	#2	-3	0	0		
		#3	0	0	0		
		#4	3	0	0		
		#5	6	0	0		
		#6	9	0	0		
		#7	12	0	0		
		#8	15	0	0		
		#9	18	0	0		
		#10	21	0	0		
		#11	24	0	0		
		#12	27	0	0		
		#13	30	0	0		
		#14	33	0	0		
		#15	36	0	0		
		#16	39	0	0		
		#17	42	0	0		
		#18	45	0	0		
		#19	48	2	2	100	
		#20	51	32	10	31.25	
		#21	54	104	46	44.23076923	
		#22	57	160	84	52.5	
		#23	60	125	71	56.8	
		#24	63	61	37	60.6557377	

F	G	H	I
107.15	84.83	85.90	
106.14	84.84	85.90	
107.12	85.86	84.83	
107.12	85.86	84.84	
107.23	85.81	-6.00	
108.00	-6.00	85.87	
108.00	86.27	85.81	
108.22	86.90	86.27	
108.64	86.26	86.90	

6. When performing MSNA analysis for CPT or EPR (or other stressors), using the DATAPAD_BEAT macro may result in incorrect data exportation in individuals with very high heart rates. If erroneous BP data is noticed, then simply open the DATAPAD_BEAT macro, and shorten the RegionSecs time period from 0.25 to 0.15, or less, depending on the HR. This may require a bit of thinking and troubleshooting, and is very rare to encounter, but something to be aware of. This error may appear as repeated BP values across several cardiac cycles.

Compiling data using R studio script. Skip this part if choosing to not use the script to compile the data.

1. First, ensure that all Excel workbooks are in one folder under one file pathway.
2. Download R and R studio.
3. The only portions that must be changed are specified below.

```
###
#copy pathway into this section of code and change all backslashes to forward slashes
#anytime you tell R to look for a file it will look in this folder
setwd("C:/Users/Miguel Anselmo/OneDrive/Desktop/PhD/Data Analysis/LabChart MSNA Analysis/REVEAL")
###
```

4.

Copy and paste the folder file pathway where the data exists. Next, **change all backslashes to forward slashes.**

```
15 #repeat this step for all participants
16 ##first specify file name then which sheet is you want the data to extracted from
17 #repeat this step for all participants
18 ##first specify file name then which sheet is you want the data to extracted from
19
20 participant1 <- read_excel("#Participant Sheet Name", sheet = "CumulativeSheet")
21 participant2 <- read_excel("#Participant Sheet Name", sheet = "CumulativeSheet")
```

5.

Next, change the Participant Sheet Name section to each participant sheet that is available. Use the exact name that is specified in the Excel Folder.

6.

```
68 cumulativeworkbook <- bind_rows(  
69   participant1, participant2,  
70   participant3, participant4,  
71   participant5, participant6,  
72   participant7, participant8,  
73   participant9, participant10,  
74   participant11, participant12,  
75   participant13, participant14,  
76   participant15, participant16,  
77   participant17, participant18,  
78   participant19, participant20,  
79   participant21, participant22,  
80   participant23, participant24,  
81   participant25, participant26,  
82   participant27, participant28,  
83   participant29, participant30,  
84   participant31, participant32  
85 )
```

Next, bind all participant files. This is an example of a dataset where there are 32 participant profiles. This and the previous step require some work, but there are methods of using AI sources to compile this code quickly.

7.

```
88 write.csv(cumulativeworkbook, "REVEAL_Cumulative_MSNA_Workbook.csv")  
89  
90 #highlight entire script before you run
```

Update the desired workbook name. Then, highlight the entire script before pressing “run.” The compiled workbook should now live in the same folder pathway where all other Excel files are housed.